# WELLCARE

## Electric Foot Warmer

Warmth in all directions | Carefully designed

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B809.101-J1024

**WELLCARE CO., LTD.**

16F, No. 738, Zhong Zheng Rd., Zhong He District, New Taipei City 235, Taiwan

Tel. 886 2 8226 2899 Fax. 886 2 8226 2889

[www.wellcare-global.com](http://www.wellcare-global.com)

## Explanation of symbols

Read instructions

Do not use when folded or creased

Do not pierce with needles

Not to be used by very young children  
(0-3 years)

Do not dry clean

Wash by hand.  
Maximum temperature 30°C

The lining can be washed in extra gentle  
washing program at 30°C!

Do not bleach

Do not tumble dry

Do not iron

**IPX1**

Drip-proof! (applicable for the foot warmer)

Thank you for purchasing the Electric Foot Warmer of Wellcare.

To enjoy the most comfortable experience and the greatest satisfaction, please take a few minutes to read this instruction carefully before using this product.




---

## **IMPORTANT SAFETY INSTRUCTIONS RETAIN FOR FUTURE USE**

### **KEEP CAREFULLY IN CASE IT IS NEEDED!**

This foot warmer is designed to be used for many years if it is properly used. Read the instruction carefully to prevent damage or danger caused by incorrect use. Keep this instruction manual at an easy-to-access location so you can refer back when necessary.

If you let somebody else use this product, please pass on this instruction manual as well and make sure the user has read it before using.
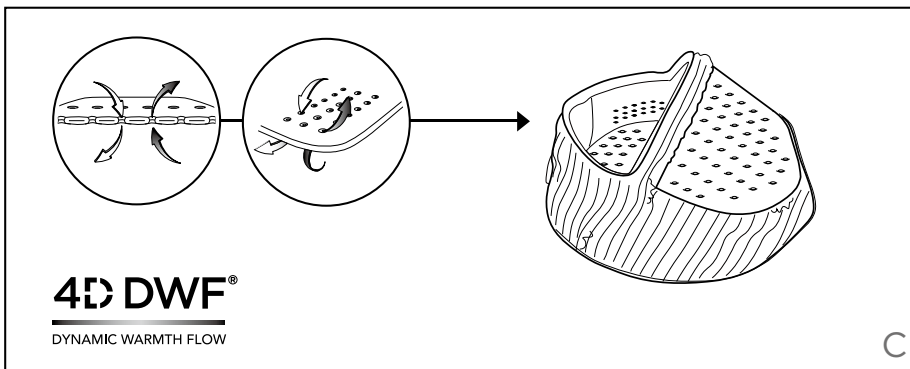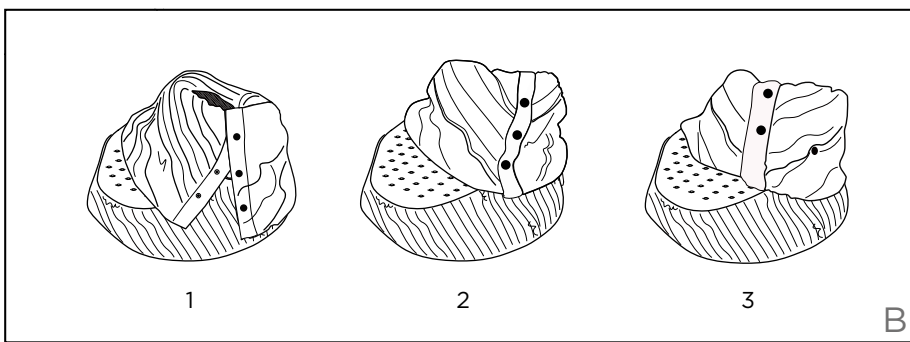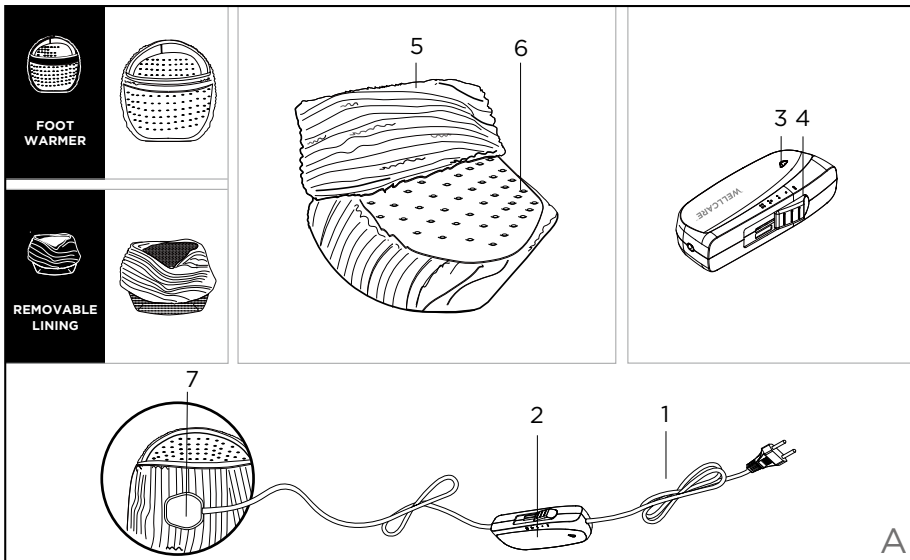
## **4D DWF® Technology** (please see fig. C)

This electric foot warmer utilizes WELLCARE's recent and the most innovative technology - 4D DWF®. While this technology pioneers in the flexible heating industry, its production process is also the state of the art

Our physical living environment consists of three dimensions, width, length, and depth; hence, the term "3D". As Wellcare introduces to you the 4D DWF® technology, the "temperature" brought to you by dynamic flow of warmth is the definition of the fourth dimension. 4D DWF® employs concepts of breathing system which creates full-body heat circulation for heat to disperse evenly, thoroughly, and quickly. The heat is generated by heating elements, builds up and gathers around several hundred uniquely designed circular holes. Through these holes, the heat flows bi-directionally to create a gentle heat flow and the thermal energy can flow rapidly to warm up human bodies.

There are many benefits of this innovative technology :

- Better breathability
- Rapid heating
- Higher temperature
- Uniquely designed holes to provide 360° warmth circulation
- Greater heat distributing area
- Effective heat transmission
- Strands of warmth flow retain air moisture
- Warmth kept at a constant temperature
- Energy saving

## **This appliance is not intended for medical use in hospitals.**

- Please read carefully and follow the instructions.
- Examine the appliance frequently for signs of wear or damage. If there are such signs or if the appliance has been misused or does not work, return it to the supplier before switching it on again.
- Do not use this appliance if there is visible damage.
- Use the appliance only as described in this instruction manual.
- Before connecting the foot warmer to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- The foot warmer must not be used when folded or crumpled.
- Do not use the appliance for persons who cannot move without help, small children, persons who cannot feel heat, or animals. Persons who can feel heat to some extent may only use the appliance with supervision.
- This appliance must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.
- Do not pierce the foot warmer with needles, safety pins or other sharp objects.
- Do not tug, twist or make sharp bends to the

cables.

- Do not carry or unplug the foot warmer by pulling its mains cable.
- Do not expose the controller, the connection box and the mains cable to moisture or liquids.
- Do not use this appliance when wet and only use it in dry environments.
- Do not operate the foot warmer unsupervised.
- Caution! Never fall asleep while the appliance is in operation.
- Prolonged use can cause burns to the skin.
- When the appliance will be used for an extended period (falling asleep), the temperature setting should be adjusted to the lowest level (setting 1) for continuous use to avoid overheating on the same body part, which may lead to burns to the skin.
- A damaged mains cable cannot be replaced. If the mains cable is damaged, the whole foot warmer must be disposed of.
- While the foot warmer is switched on,
  - Do not place any large or heavy objects on it such as suitcases, etc.,
  - Do not place any heat sources on it such as hot water bottles, heating pads or such like.
  - Due to the components in the controller will make it slightly warm, never let the controller to be covered or placed on the heating pad, pillows and the like or lie on the electric mattress heater.

- Do not fold the foot warmer when it is hot.
- Children should be supervised to ensure that they do not play with this appliance.
- Children below the age of three (0-3 years old) are not allowed to use this appliance because they are not able to react to overheating.
- This appliance cannot be used by young children over the age of three (3 ~ 8 years old) unless the controller has been preset by a parent or guardian or the child has been instructed and how to operate the blanket safely.
- This appliance may be used by children over the age of eight (8), and persons with reduced physical, sensory or mental capabilities or lack of experience or knowledge if they have been supervised or instructed concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be made by children without supervision.
- If you have any concerns about health issues, please consult with your doctor before using the appliance.
- Do not use the foot warmer on parts of the body that are swollen, injured or inflamed.
- If you experience any pain or discomfort while using the unit, stop using it immediately.
- If you use a pacemaker, please consult your doctor and the manufacturer of the pacemaker before using this product due to the EMF

(electromagnetic fields) emanated from this electric product may interfere with functions of your pacemaker under certain circumstances.

- This foot warmer is only designed to warm up the human body. Do not use it on pets or other animals.
- When the foot warmer is not in use:
  - Set the slide switch on the controller to setting “0” (OFF)
  - The plug should be removed from the wall socket and the foot warmer should be stored in a dry place, without additional weight on it.
- When the foot warmer will not be used for a period of time, we recommend that you store it in the original packaging.
  - Allow it to cool down before the storage.
  - Do not crease the foot warmer by placing items on top of it during storage.
  - Store it in a dry place away from any liquids or sources of heat.
- Do not wear shoes when using the foot warmer because this will lower your heat-sensitivity and will make the foot warmer dirty.

## **ITEMS INCLUDED IN THE PACKAGE**

- 1 Foot Warmer
- 1 Removable Fleece Lining
- 1 Controller with Cable and Power Plug
- 1 Instruction Manual
- 1 Product Catalogue
- 1 Service Instruction
- 1 Membership Card

## **DESCRIPTION OF THE DEVICE AND CONTROLLER** (please see fig.A)

1. Mains cable with plug
2. Controller
3. APS (Active Protection System) with 3-colored LED display (on both sides)
4. Slide switch for ON, OFF and temperature settings (1, 2, 3, 4)
5. Removable fleece lining
6. Foot warmer
7. Connection box

## **ADJUSTABLE FLEECE LINING** (please see fig.B)

The snap button design on the fleece lining can help you to adjust the lining to wrap your feet higher, cover more area and feel snuggler.

## **PREPARING FOR USE**

1. Always make sure the foot warmer (6) is

without any folds, creases or ruck-up.

2. Always use the foot warmer (6) with the removable fleece lining (5) and only use the one supplied with it.
3. Make sure the controller (2) is not covered by the sheet, pillow or any objects.
4. Plug the mains cable (1) to the wall socket.  
Unplug the appliance from the wall socket when it is no longer in use.

## SWITCHING ON

---

Before each time you use the foot warmer (6), please make sure it is not bundled, folded or creased before it is switched on.

Connect the product only to 230VAC. Switch on the product by using the hand-operated slide switch (4) and the light on the APS (3) will be turned on to ensure the product has been operated correctly.

### **Temperature settings :**

**4 = maximum warmth**

**2-3 = moderate warmth**

**1 = minimum warmth**

**0 = OFF**

In order to warm up the foot warmer faster, it is recommended to initially set the foot warmer at the highest temperature setting.

## **APS (Active Protection System) – 3-Color LED**

This product is equipped with APS (Active Protection System) (3) technology which detects errors on the operation or component failures. By the 3-color LEDs which can be seen from both sides of the controller (2), APS can help you to operate the foot warmer easily and safely.

Green – The product has reached the desired temperature level.

Orange – The product is heating up to reach the desired temperature level.

Due to the temperature variation in different environments (room, bed, etc.), APS LED may change between green and orange colors within seconds which indicates continuous heating when the product is in use. Under this circumstance, please note the product is functioning properly and normally.

Red light – Possible malfunction or overheating, switch off the appliance immediately.

In the following situations, the red light will come on to indicate warning:


- When the critical components are damaged: Please disconnect the mains cable from wall outlet immediately.
- When the heating wires or other components are broken and the device cannot work anymore.

The warnings mentioned above cannot cover all defective situations as the damages on some

components will not lead to red light.


### TURBO HEATING

The foot warmer is with a turbo heating function for user to feel warmth within 7-10 minutes after it is switched on.
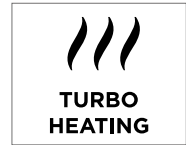
### TIMER FUNCTION

Now you can enjoy the pleasant warmth especially for your feet. This foot warmer is equipped with an automatic switch-off function. After approx. 90 minutes, the foot warmer will be switched off automatically.

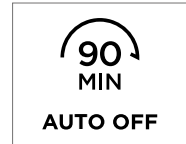In order to put the foot warmer into operation again you must set the slide switch (4) to setting '0' first.

After about 5 seconds, you can switch on the foot warmer again.

### CAUTION

- Do not use the foot warmer when it is folded, creased or drawn together.
- The electronic components inside the controller may cause the controller slightly warm when the product is in use. Therefore, the controller must not be covered or placed on top of any heating sources when it is in operation.
- You should not gauge the foot warmer's temperature by hands because hands are less sensitive to heat than the skin on the rest of the body parts. The temperature may feel lukewarm by hands but could be unpleasantly hot to feet.

### CLEANING AND CARE

- Before cleaning the foot warmer (6), always unplug the product from the wall socket first. Otherwise, there is a risk of electric shock.
- To clean the controller (2), use a dry, lint-free cloth. Do not use any chemical or abrasive cleaning detergents.
- Please note that the foot warmer (6) should not be dry-cleaned, wrung out, machine-dried, pressed, ironed or bleached.
- The foot warmer (6) can be hand washed at 30°C; small marks on the foot warmer can be cleaned with a slightly damp cloth or a soft brush. The connection box (7) must not come into contact with water or any liquids.
- The removable fleece lining (5) is machine washable. Set washing

machine to a gentle washing cycle at 30°C, the fleece lining should only be washed together with other textiles. Use a mild washing detergent for delicate fabrics and dose it according to the manufacturer's instructions. Always use a laundry bag for better protection of the fabric.

- For dark color products, always wash separately.
- After washing, reshape the foot warmer (6) and the fleece lining (5) to its original dimensions while it is still damp, and lay them flat on a moisture-absorbent towel, mat or clothes horse to dry it.
- Do not use pegs or similar items to attach the foot warmer (6) to the clothes horse, they may damage the foot warmer.
- Only use the foot warmer (6) and the fleece lining (5) when they are complete dry.
- Never switch on the foot warmer to dry it.

## STORAGE

If the foot warmer will not be used for a long time, we highly recommend:

- Always let the appliance to cool down first before the storage.
- Store it in the original packaging and put it in a dry place.
- Do not put any object on the top of the appliance.

## DISPOSAL

This appliance must not be disposed together with domestic waste but at a suitable local collection or recycling point. Please dispose of the appliance in accordance with Directive on Waste Electrical and Electronic Equipment (WEEE). If there is any question, please contact the municipal authority for more information about disposal in your area.

## WARRANTY

Please note Wellcare provides two-year limited warranty for this product. Wellcare warrants that this product will be free from defects in materials and workmanship for a period of two years beginning on the date of purchase.

During the two-year period beginning on the date of purchase, Wellcare will replace defective products covered under this limited warranty. Please note a clear copy of the purchase receipt for the product is required. To request limited warranty service, please return the defective

product to the store where you purchased it with the clear copy of purchase receipt within the limited warranty period. The replacement product will assume the remaining limited warranty of the original product.

We advise you to call our helpdesk for any problems regarding the product before returning it to the store where you purchased it.

However, please note that our limited warranty does not cover any damage due to shipment of the products to you and external causes, including accident, abuse, misuse, problems with electrical power, usage not in accordance with product instructions, failure to perform reasonable maintenance, problems caused by use of parts and components not supplied by Wellcare, and problems caused by unauthorized repair of Wellcare's product. Furthermore, our limited warranty does not cover any damage due to the defects known to exist by you when you purchase the product.

This limited warranty gives you specific legal rights, and you may also have other rights that vary from country to country (or jurisdiction to jurisdiction). Wellcare's responsibility for malfunctions and defects in product is limited to replacement as set forth in this limited warranty statement. All express and implied warranties for the product, including but not limited to any implied warranties of an conditions of merchantability and fitness for a particular purpose, are limited in duration to the limited warranty period set forth above and no warranties, whether express or implied, will apply after such period. Wellcare does not accept liability beyond the remedies set forth in this limited warranty statement or liability for indirect, incidental or consequential damages.

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